**MANUALE D'USO**  
**USER'S MANUAL**

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**AUTONOMOUS POSITIONER MYBACK**

Codice/Code	3957_2021
Descrizione/ Description	Myback Eco -Fit
Lingua/Language	Italiano / Inglese
Redatto da/ Written by	Donati Srl
Rev.	1
Data /Date	01-04-2021

# SOMMARIO / SUMMARY

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----

Thanks for having chosen Myback ECO FIT.

This manual contains some information concerning the technical performances of the equipment, its functioning, its maintenance and its safety use.

The manual is addressed to all the users: at any level they (customers, operators, technicians) are obliged to read carefully the manual before doing any operation on the equipment.

The manual must be kept in proximity of the equipment, available immediately to all the users and protected from any possible damage which could affect it.

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----

Myback's purpose is to allow the autonomous implementation of a stretching and back discharge exercise through the gentle lifting of the pelvis with back and legs outstretched

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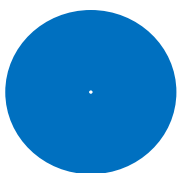
**AFTER USE ALL WASTE ELECTRICAL AND ELECTRONIC EQUIPMENT MUST  
BE TRANSPORTED TO AUTHORISED TREATMENT FACILITIES.**

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### **Myback ECO FIT Code 3957\_2021**

complies with the essential requirements about safety and health provided by machineries directive 2006/42/CE. It complies with the UNI EN ISO 20957-1 2013 standard (“Fixed device for training, Part 1: (general requirements and testing methods)

CLASS OF USE: H

MAXIMUM USER WEIGHT: 130 kg

It also declares that :

The electric movement device, complying with the standards :

EN 60635-1

The technical documentation issued by the producer is kept in his company

<b>CE</b>	<b>Donati Srl</b> Via Francesca Nord, 50 56010 Vicopisano (Pi) Italy Tel.+39 050 796022	
Peso max utilizzatore/ Max user weight 120 Kg Class H Norme UNI ISO 20957-1 2013		
Mod.	3957_2021	
Anno/year	xxxx	
Lotto/Lot	xxx	

*Timbro e firma rappresentante Legale/  
Signature of legal representative*

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Declaration of conformity of the electrical components inside the equipment

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Shiyong Minyin Industrial Zone, Hengli Town, Dongguan City, Guangdong, China.

### EU Declaration of Conformity

According to Directive 2014/30/EU

Hereby we declare that the concept and design of the following listed products, as well as the type design supplied by Timotion, conform to basic requirements of the above mentioned EMC directive. And modification or change of the product without our approval will void the validity of this declaration.

Manufacturer: Timotion Technology Co. Ltd.  
Shiyong Minyin Industrial Zone, Hengli Town,  
Dongguan City, Guangdong, China

From date of construction: 4 Nov. 2016

#### Products:

Actuators: TA Series, MA Series, JP Series, VN Series.

Control Boxes: TC Series, TBB Series.

Handsets: TH Series, TMH Series.

Columns: TL Series, TVL Series

Gear Motor: TGM Series.

Foot Switches: TFS Series.

Cables: TYC, TEC.

Battery Charger: TBC Series.

Bottom Light: TBL.

Safety Strip: TSS.

RF Receiver: TRF Series.

Adjustable Frame: TEK Series

Power Supplies: TP Series.

Nursing Pad: TNP Series.

Junction Box: TJB Series.

Massage Motors: TMM Series, TAM Series.

Furniture Handsets: TFH Series, TDH Series.

Heating Pad: THP.

Lock Boxes: TLB.

External Switch: TES.

Reading Light: TRL.

Emergency Button: TEB Series

Other: TNC, 3\*9V Battery Box, TFA Series

SP Series, TCS1, TWD Series, TWC Series

#### Applied Harmonized standards:

EN 55014-1:2017/A11:2020, EN 55014-2:2015, EN 55011:2016+A1:2017+A11:2020

EN IEC 61000-3-2:2019, EN 61000-3-3:2013+A1:2019, EN IEC 61000-6-1:2019,

EN 61000-6-3: 2007+A1:2011

EN IEC 61000-6-4:2019, EN IEC 61000-6-2:2019

IEC 60601-1-2:2014, EN 60601-1-2:2015

Jason Lin / Deputy Managing Director

## EU Declaration of Conformity

According to Directive 2014/35/EU Low Voltage Directive

Hereby we declare that the concept and design of the following listed products, as well as the type c design supplied by Timotion, conform to basic requirements of the above mentioned LVD directives. An modification or change of the product without our approval will void the validity of this declaration.

Manufacturer: Timotion Technology Co. Ltd.  
Shiyong Minying Industrial Zone, Hengli Town,  
Dongguan City, GuangDong, China


From date of construction: 20 Apr. 2016

Products:

Actuators: TA Series, MA Series, JP Series, VN Series.	Power Supplies: TP Series.
Control Boxes: TC Series, TBB Series.	Nursing Pad: TNP Series.
Handsets: TH Series, TMH Series.	Junction Box: TJB Series.
Columns: TL Series. TVL Series	Massage Motors: TMM Series, TAM Series.
Gear Motor: TGM Series.	Furniture Handsets: TFH Series, TDH Series.
Foot Switches: TFS Series.	Heating Pad: THP.
Cables: TYC, TEC.	Lock Boxes: TLB.
Battery Charger: TBC Series.	External Switch: TES.
Bottom Light: TBL.	Reading Light: TRL.
Safety Strip: TSS.	Emergency Button: TEB Series
RF Receiver: TRF Series.	Other: TNC, 3*9V Battery Box, TFA Series
Adjustable Frame: TEK Series	SP Series, TCS1, TWD Series, TWC Series

Applied Harmonized standards:

EN60335-1:A1:2019+A2:2019+A14:2019, EN60335-2-17:2013/A1:2020,  
EN 60335-2-103:2015, EN 62233:2008, EN/IEC 61558-1:2019  
EN 61558-2-16:2009+A1:2013  
IEC 60601-1:2005+A1:2012 (3rd Edition)

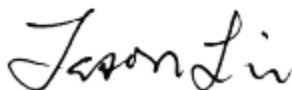Jason Lin / Deputy Managing Director

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
----

Myback ECO FIT consists of a vertically moving column, electrically operated by a remote control, and supports for legs and ankles:

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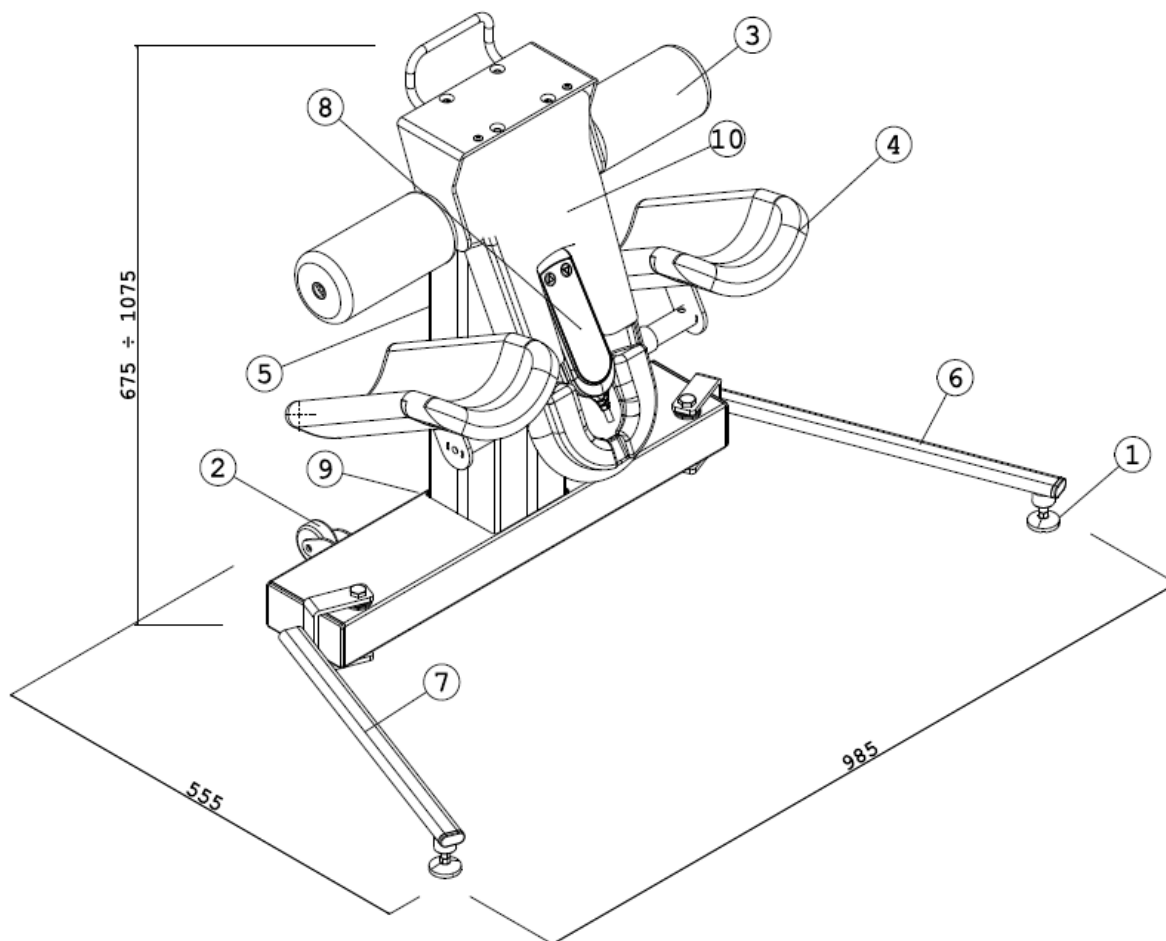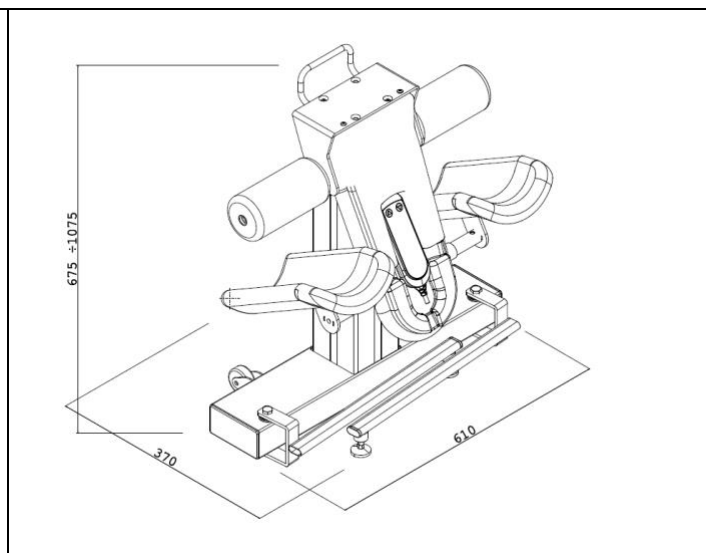
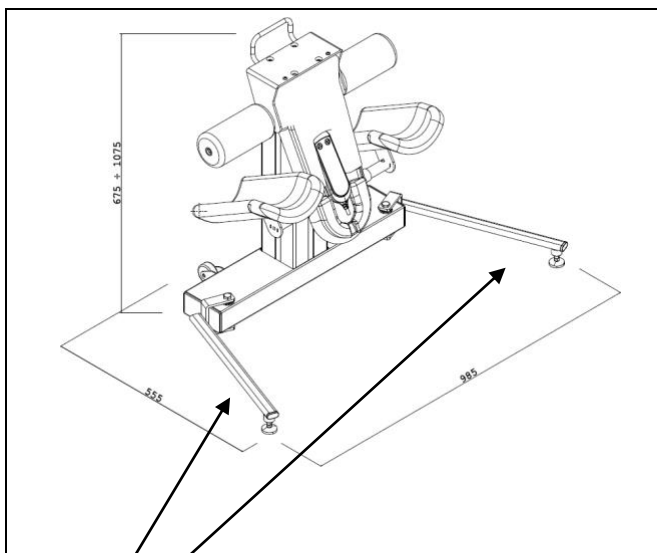
**It's forbidden to do the exercise with the supports in the closed position**

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Before using the equipment, please read this user's manual carefully in all its parts, in particular the warning instructions specified in paragraph 2.

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The autonomous positioner must be used by a single person who enjoys full physical capabilities and uses it in the correct way, as shown in the picture.

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The exercise should be performed with the supports rotated forward as shown (Fig B), otherwise the equipment could turn over

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The equipment is not a toy and cannot be used by children under 16 years old unless under the direct control of a paediatrician.

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Do not leave a child or a disabled person unattended on Myback ECO FIT equipment.

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Myback ECO FIT should use only for its intended use, that is as a positioner to assume a position of lumbar relaxation; it's not a workout bench and should not be used for endurance or speed training. The standard workout is: 2 minutes of exercises and 5 min of rest.

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Do not connect the power cable to nets having different voltage than those indicated in the present manual .

Do not wet the positioner (it contains electrical parts)

Do not dismount or tamper with the components, do not remove the fences.

Do not pull, cut or touch with acuminated objects the connecting cables.

Do not expose to direct sunrays the remote control.

Do not use the equipment or its components to different uses from those indicated in this manual.

Do not use My Back if damaged, overheated or in presence of smoke

The myback is equipped with an auxiliary battery useful for lowering the column in the event of a power failure

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The equipment is delivered fully assembled, the packaging of Myback allows to minimize the risks coming from the conditions of transport. It also has handles that allow you to carry it to a desired place, provided the conditions hereunder described

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Handling of Myback when not packed:

Instructions for its use through the top handle:

- Ensure that the swing bases are in the closed position as shown in paragraph
- Stand at the back of the equipment and hold the top handle.
- Holding the handle incline the equipment (towards your direction) about 30° until you notice that the back swing bases start to rise from the ground and use the wheels (Fig A pos.5).
- Keeping the equipment inclined, you can drag it or short times on the same plane.
- Do not drop rudely the equipment while folding

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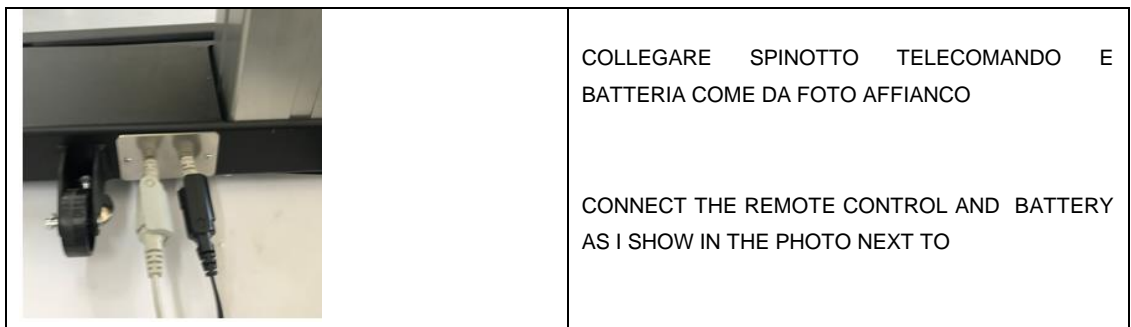
Position your MyBack on a flat level surface, in a closed environment where you have enough space (at least 1.5 meters between the equipment and the surrounding objects), to perform the exercises (lumbar stretching).

Myback must be placed near a 220V electrical socket. Before the electrical connection, verify that the supply voltage on the label is the same of the country where Myback is used.

Be sure that the connection cable does not interfere with the use of the equipment. The electrical cables must be protected from damages and humidity.

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Myback has been designed, produced and tested for functioning indoors, with adequate light, change of air, low percentage of humidity and low quantity of dust.

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Once verified that the content is as specified in this manual, proceed as follows:

- Remove the protective packaging materials used for the transport, ensure that the conditions of the machine are suitable for the proper functioning of the same, then check that all parts of the machine are in good conditions before its use.
- Keep the packaging material for the validation of warranty.
- Lift the tool from the box in 2 people to avoid excessive strain on the back

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Remote control: On t e control panel only two commands are available. Push the button  for sliding up the support . Push the button  for sliding down the support

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Place a fitness mat near the front of the equipment  
(the mat is not included with the Myback)

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
THE LATERAL SUPPORTS MUST BE IN OPEN POSITIONS (SEE PICTURE. B)  
IT IS STRICTLY FORBIDDEN TO USE MYBACK WITH FOLDED (ie in its position when paked)

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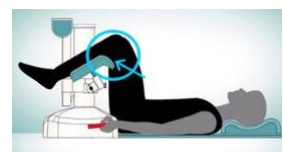Sit on the mat by placing your head on a thin cushion, with legs on either side of the column.

Place your legs on the two lower supports, with your popliteal fossa (the back of the knee joint) in the correct position, as shown on the right.

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Pushing the button “▲” on the remote control you can raise your pelvis to the position of "discharge" you want. The rear cylinders hold the ankles and prevent the legs from slipping back. The final height of lift varies according to user's characteristics, generally it's recommended to lift the pelvis from the mat of 4-5 cm (the height of a punch as an indication). Once reached the position of passive discharge, keep it for an indicative time of 2-10 minutes, or until you feel the slightest discomfort or fatigue. (see the advice section for the exercises performed in a position of discharge)

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After the expected time, lower completely the full support through the button “▼” on the remote control.

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Remove the legs from the support and shifting ankles and feet support towards the outer side, then place your feet on the ground as in the starting position.

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In some cases of longer time of positioning in lumbar discharge, you could feel a slight pins and needles on the popliteal cavity and on the feet. This sensation disappears few minutes later after having concluded the session. It is better to stop the session and start again after some minutes.

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BEFORE GETTING UP put the remote control on the column's central support

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

In case of power outage while using Myback and the auxiliary battery was discharged , it will not be possible to lift or lower the leg's support.

Follow the instructions below to get out of the positioner :

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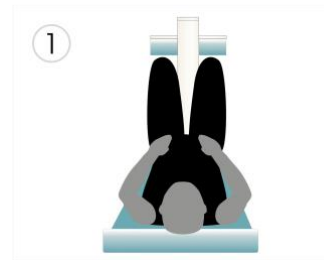From the raised position (Picture 1), free your ankles by moving one foot at a time towards the outside (Picture. 2), and slide back placing your feet as in the starting position

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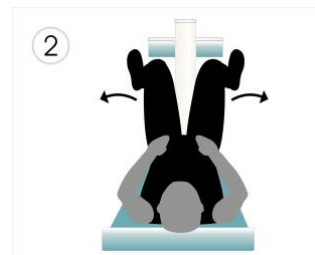### Removal of the emergency auxiliary batteries

Operation to be done with the tool disconnected from the mains

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Remove the screw using a screwdriver, cross head

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Replace the three batteries 9 volt and re-insert in the order as per original installation

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MyBack does not require a periodic maintenance;

All the cleaning and maintenance operations must be done when the equipment is not functioning

Please note that only qualified technicians can work on the electrical parts.

Clean the exterior of MyBack using a dry cloth or slightly moistened with mild detergent, periodically clean the column with a dry cloth or a brush.

In the event that the power or connection cables to the panel are damaged, they must be replaced by the manufacturer or by his technical support

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MyBack ECO\_FIT is guaranteed for a period of 24 months from the date of delivery or proof of purchase.

Donati s.r.l. will repair or replace free of charge or at its own judgment all products or components that show manufacturing defects or malfunction, provided that they are reported in writing within the warranty period.

The warranty is void for those products tampered with or damaged by causes not attributable to Donati Srl, as well as in the case of use or installation not in accordance with the instructions contained in the user's manual supplied with the product.

The warranty does not cover any damage due to inadequate maintenance or repairs made without written authorization from the manufacturer.

The warranty is limited to repairing the defective product at our office, or on the premises of the Customer, or to its replacement by Donati Srl's sole opinion.

Nothing is due to the buyer for the period during which he will be left without equipment, no form of compensation shall be recognised under any circumstances

Note:

- Any other warranty is not under our responsibility.
- This warranty does not affect your legal rights
-

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---

The criteria below have been followed for the choice of manufacturing solutions :

- a) removal or reduction of the risks in the best possible way, integrating the security problems during the drawing of the project and the manufacturing of the equipment.
- b) adoption of appropriate measures of protection where the risks could not be eliminated.
- c) Information to the users of residual risks due to possible ineffectiveness of the protection cautions adopted.

All has been considered not only in normal functioning conditions but also in abnormal functioning conditions even if possible .

The equipment has been projected and produced following the ergonomics principals for reducing as much as possible the troubles, the strain and the stress of the user .

The user's manual contains the instructions for handling. The materials used for the construction of the equipment are not dangerous for the security and the health of the people.

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The control devices are projected and made for being safe and reliable.

The electrical control devices are visible, immediately accessible to the user and marked with indications for easy interpretation.

Power outages and their restoration do not cause any dangerous situation.

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The instructions contain indications to prevent an inappropriate use of the equipment.

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All the components of the equipment are designed and manufactured so that, in the intended operating conditions, they are stable for being used without any risk of overturning and falling.

The elements have sufficient strength and proper use of the environment provided by the manufacturer.

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Fixed protections are used for all parts of the equipment presenting hazards

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The equipment has been designed, manufactured and equipped to prevent the risks connected with electrical power, according to the actual laws .

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No flammable materials have been used.

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The equipment works under the noise limit of 70 dB(A)

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The equipment does not produce any vibration

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Chemical essences and compounds are not requested for normal daily cleaning of the equipment.

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Maintenance, lubrication, cleaning and repairing operations can be done when the equipment is not functioning.

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All the necessary information for the use of the equipment and its risks are clearly reported on the user's manual .

Devices for personal protection are not provided

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Myback is an electric equipment and must be correctly disposed. In particular, its collection as waste must be done by a specialised company duly authorized to the transportation and disposal of electric and electronic devices.

### INFORMATIVE REPORT TO THE USERS OF AEE

According to art.13 concerning the decree law dated July 25 2005 no. 151 "Application of Directives 2002/95/CE,2002/96/CE and 2003/108/CE, regarding the reduction of dangerous substances in electric and electronic equipment as well as waste disposal".

This product is conform to the directive EU 2002/96/EC.

The components of the product , marked with the crossed symbol of the container for waste reported on the equipment, are subjected to the compulsory waste selection at the end of their cycle, in all the countries belonging to the European Community and in those who adopt the separate waste collection :

- It is forbidden to dispose such components as urban undifferentiated wastes and cannot be disposed together other domestic wastes.
- It is compulsory to collect separately such components and give them to the authorised treatment facilities or to return them to an authorized distributor when buying an equivalent product.

#### ATTENTION !!!!

- An unauthorized disposal of the product may cause damages to the environment and involves the application of administrative sanctions foreseen by the existing set of rules

- It is forbidden to use the electric and electronic components of MyBack® in inappropriate ways, it is forbidden to dismount and use them if damaged. Such components could have harmful health effects.

What above described is referred exclusively to all the components of Myback marked with the symbol

For additional information about the correct disposal of the product and the available waste collection systems, please contact your local treatment facility or contact Donati s.r.l. at +39 050 796022 , e.mail [info@myback.it](mailto:info@myback.it)

**National Register of Producers AEE**